

PORODO
LIFESTYLE



Porodo Lifestyle

**Lavender
Ice Cold IPL Hair Removal**

SKU: PDLFSTY11WH

Table of Contents

Product Overview	2
Specifications	2
Safety Precautions	3
Schematic View	3
Package Contents	4
Device Operation	5
Safety Guidelines	6
Operation & Caution Guidelines	8
Restrictions on Use	11
Skin Sensitivity Test	11
Light Intensity Settings	11
Pre-Use Preparation	12
Instructions for Use	13
Treatment Schedule Guide	14
Post-Treatment Care	14
Cleaning & Maintenance	15
FAQ	15
Disposal	17
Warranty	18
Contact Us	18

Product Overview

The Porodo Lavender Ice Cold IPL Hair Removal device offers effective hair removal with 1-9 adjustable intensity levels and multiple flash modes for Hair Removal (HR), Skin Care (SC), and Regeneration (RA). It features a cooling function for added comfort during treatment and includes automatic flash and ice cooling modes for ease of use. This compact device provides professional-quality results for smooth, flawless skin at home.


Specifications


Rated Voltage	100V-240V
Rated Frequency	50/60Hz
Power Cord	UK 3-Pin Plug
Rated Power	40W
Light Energy	10-15J/cm2
Wavelength Range	470nm to 1200nm
Energy levels	9 Adjustable Modes
Operation Modes	Hair Removal, Skin Care, Rejuvenation/Anti-Aging
Storage Temperature	0°C to 40°C
Operating Temperature	5°C to 35°C
Product Weight	358g
Model Number	PDLFSTY11

Safety Precautions


1. Please use the product as instructed in the manual to ensure proper use and safety.
2. This device is not for use by minors under the age of 16. Keep it out of reach of children to avoid any safety risks.
3. In line with our commitment to sustainable development, we reserve the right to update or modify our products, including changes to their function, interface, or appearance.
4. Due to product upgrades, if there is any difference between the description and the actual product, the details of the real product will apply.
5. For your safety, always wear protective glasses when using this device and keep the flash window away from your eyes during use.

Schematic View

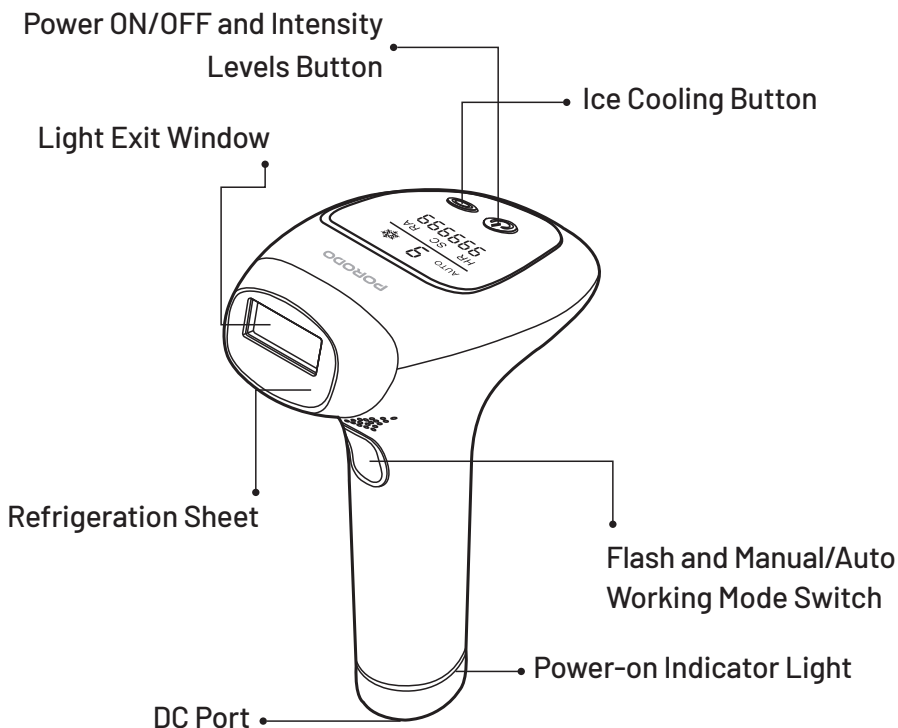


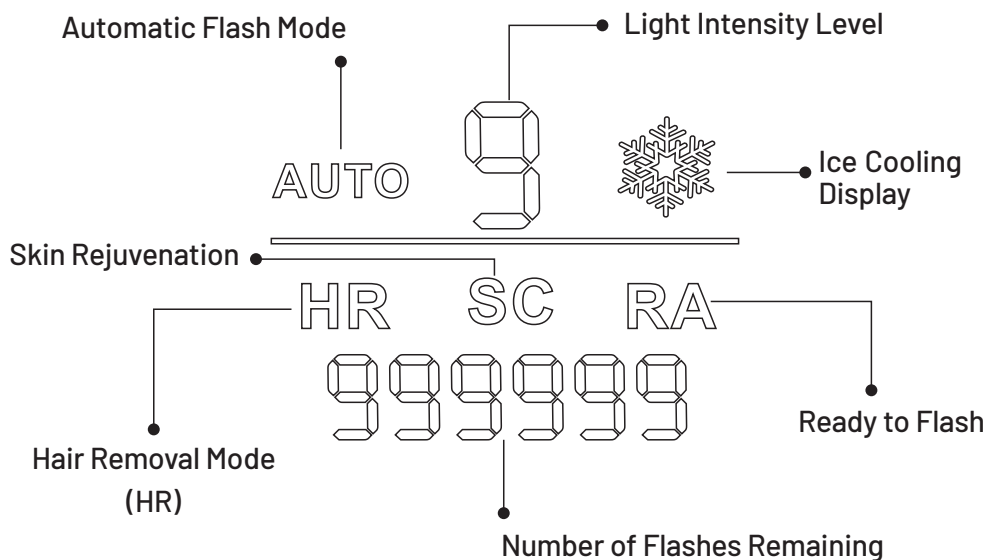


Package Contents

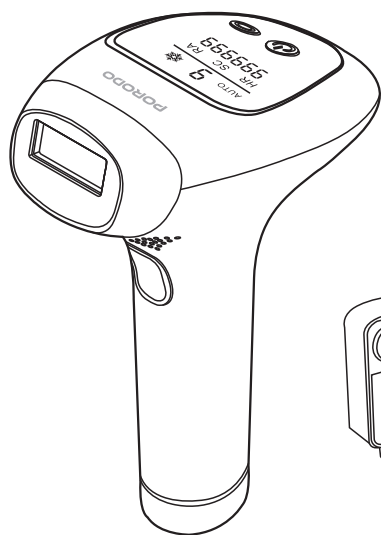

Product (Hair Removal Device) ×1

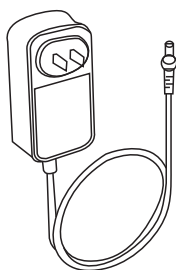

The Power Adapter ×1

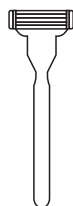

Razor ×1



Goggles ×1

Device Operation

a. Power On:

- 1.** After connecting the device to a power supply, the power-on indicator light will turn on.
- 2.** To turn on the machine, press the Power button briefly or long press.
- 3.** Once the device is on, press the Intensity Levels button to adjust the intensity from levels 1 to 9.
- 4.** The screen will display the selected intensity level accordingly.
Levels 1-2: Skin Rejuvenation (SC) mode is displayed.
Levels 3-9: Hair Removal (HR) mode is displayed.

b. Power Off:

To turn off the device, long press the Power button for 3 seconds. If the device is left on without any operation for 15 minutes, it will automatically shut down for safety.

c. Flash Mode Functions

1. Manual Flash Mode:

Manual flash mode is activated automatically when the power is on.

Place the device vertically on your skin and briefly press the Flash button to begin your IPL treatment.

(Note: No flash will be emitted if the device is not in full contact with the skin.)

2. Auto Flash Mode:

After pressing the Flash button for 3 seconds, the "AUTO" icon will appear on the screen, indicating the device is in auto flash mode. It will automatically flash as soon as it makes full contact with the skin.

To switch back to Manual mode, press the Flash button for 3 seconds.

d. Ice Cooling Mode

To activate Ice Cooling mode, briefly press the Ice Cooling button, and the ice cooling icon will appear on the display screen.

The Ice Cooling mode can be used independently to soothe the skin. However, do not press the Flash button when using the cooling function alone to prevent skin damage.

The Manual Flash mode or Auto Flash mode can be used together with the Ice

Cooling mode. (Note: No flash will be emitted if the device is not in full contact with the skin, but the cooling function will still be active.)

e. Ice Cooling Mode + Manual Flash Mode:

The device should be placed vertically on the skin. The flash will emit with each press when the light outlet is fully attached to the skin, and the cooling effect will be activated simultaneously to relieve any skin discomfort.

f. Ice Cooling Mode + Auto Flash Mode:

The "AUTO" icon will appear on the screen after pressing the Flash button for 3 seconds, indicating the device is in Auto Flash mode. The flash will emit automatically once the light outlet is in full contact with the skin, and the cooling effect will work simultaneously to soothe the skin.

It is normal for water droplets to appear on the cooling area of the lamp when in Ice Cooling mode. After turning off the device, use a soft cloth to wipe off the water droplets to avoid false flashing.

Safety Guidelines

- 1. Usage Restrictions:** The use of this product may cause discomfort or adverse effects in individuals with allergies, skin diseases, allergic dermatitis, or sensitive skin. If any adverse effects occur during use, stop using the product immediately.
- 2.** Do not use this product on areas with compromised skin integrity, such as burns, blisters, cuts, wounds, or areas with active skin disease or infection.
- 3.** If you are taking photosensitive drugs for medical treatments, it is important to consult a doctor before using the product. Using the product without consultation may damage the skin.
- 4.** This product is not suitable for pregnant women, lactating women, individuals experiencing menstruation, those with body scars, diabetes, fragile skin, or light-sensitive skin.
- 5.** Individuals with severe heart disease, moderate to severe hypertension, malignant tumors, folliculitis, infectious skin diseases, burned skin, or anticoagulant disorders should refrain from using this product.
- 6.** Do not use if you have abnormal collagen metabolism, including a history of scars or poor wound healing.

- 7.** This product should not be used by individuals with vascular diseases, such as varicose veins or vasodilation in the treatment area.
- 8.** Avoid use if your skin is sensitive to light or prone to rashes or allergies.
- 9.** The device should not be used on areas that are infected, or by individuals with herpes simplex, skin lesions, or hematomas.
- 10.** Individuals with suppressed immune response diseases, including HIV infection, lupus erythematosus, porphyria, etc., should not use this product.
- 11.** Individuals with bleeding disorders, or those using medications such as tretinoin, fruit acid (AHA), salicylic acid (BHA), or azelaic acid, should avoid using this product.
- 12.** Use this product with caution if you have recently undergone chemical hair removal, used skin whitening products, or used vitamin A ingredients within the past 30 days.
- 13.** Do not use with essential oil products at the same time. Make sure the skin in the treatment area is clean and dry, without any lotions, lubricants, or cosmetics applied.
- 14.** If the skin is not clean and dry, it may affect the product's performance and could damage the device.
- 15.** This product should not be used by individuals who have recently received injections of fillers or toxins.
- 16.** Do not use the device on sunburned or burned skin until fully healed, as it may worsen the condition. It is also not suitable for individuals with epilepsy or those who experience light-induced seizures or migraines, as the light may make these conditions worse.
- 17.** Avoid using the device on areas with tattoos, dark spots (such as dark freckles, birthmarks, moles, or warts). These areas absorb more light, which increases the risk of skin burns.
- 18.** Do not use the device around the head, ears, nipples, genitals, or anus, as this may increase the risk of skin burns.
- 19.** If you have skin conditions such as psoriasis, eczema, herpes, wounds, or rashes, do not use the device. These conditions make the skin more sensitive and may cause swelling, blistering, or infection. Please consult a doctor to ensure the condition is treated before use.

20. Do not use this product if you have scars, light-sensitive skin, or if you are taking medications that increase photosensitivity. This includes non-steroidal anti-inflammatory agents (such as aspirin, ibuprofen, acetaminophen), tetracycline, phenothiazole, diuretics, sulfonamides, triazine imidazoleum, fluorodiazine, vinblastine, griseofulvin acid, carotene acid, all-trans vitamin A acid, isotretinoin, topical retinoic acid, or other non-steroidal anti-inflammatory agents.

21. If you have a history of light-sensitive diseases, such as porphyria, polymorphous solar rash, or solar urticaria, it is forbidden to use the device, as it may increase skin sensitivity.

Operation & Caution Guidelines

- 1.** Always wear protective glasses when operating the device.
- 2.** Before using the device, clean the areas of the body that need to be treated. Do not use any cosmetics or skincare products before use, as these may damage the device and affect its performance.
- 3.** When using the device for the first time, start with the lowest intensity level and select the appropriate setting to avoid damaging the skin. If higher intensity is used, the skin may become red or swollen. In such cases, immediately apply cooling to the skin.
- 4.** After use, ensure to replenish moisture to the skin and protect it from strong light and ultraviolet rays by applying sun protection.
- 5.** It is recommended to use pure plant-based moisturizing and repairing products for post-treatment care.
- 6.** Clean the surface of the device using a dry soft cloth or a wet cloth without alcohol.
- 7.** Do not expose the device or its accessories to liquids or humid environments.
- 8.** Keep the device in a well-ventilated, dry environment, and ensure it is kept out of reach of children.
- 9.** Do not insert metals into the cooling hole, as this may pose an electric shock hazard.

- 10.** If the product casing is damaged or accidentally filled with water, send the device to the manufacturer for repair and ensure it is safe to reuse. If water enters during operation, immediately unplug the power supply to ensure safety.
- 11.** Do not attempt to disassemble or repair the device without authorization.
- 12.** If the device malfunctions, send it back to the manufacturer with all included accessories for professional maintenance.

Warning

- 1.** Do not direct the device towards the eyes, as the strong light emitted by the product may cause damage to the eyes.
- 2.** Under normal operating conditions, avoid direct contact with the heat dissipation outlet, as the temperature of the outlet may cause burns to the skin.
- 3.** Under normal operating conditions, do not insert any objects into the device, as this may cause damage to the human body.
- 4.** Avoid applying antiperspirant immediately after treating the underarms.
- 5.** To prevent scalding or malfunction, regularly check and clean the irradiation port to remove any dirt or debris.
- 6.** Do not place the irradiation port near items that are easily demagnetized, as this may cause device failure and damage.
- 7.** The device should not be used if the vent is blocked, as this could lead to smoke or fire.
- 8.** For proper operation, avoid pointing the irradiation port at black objects.
- 9.** To prevent damage such as blackening, discoloration, or breakage, do not irradiate clothes, hair, or other items.
- 10.** The device is not intended for use on pets, as it may cause injury.
- 11.** The device may cause eye damage (including potential blindness) or skin damage if the instructions are not properly followed. Always read and adhere to the instructions provided.
- 12.** Do not use the device on surfaces that require excessive use or may result in excessive pulse stacking and/or frequency, as this could damage the skin.
- 13.** Misuse of the device may cause damage to the eyes. Always protect your eyes from exposure and avoid using the device near the eyes.
- 14.** Before using the device, perform a skin test and confirm the skin condition after 2 hours. If there is no abnormality (such as redness or blisters) in the treated area, the product can be used.

- 15.** During use, ensure to wear the goggles and place the lamp window close to the skin for the depilation process.
- 16.** After use, wipe the hair and skincare products attached to the device to avoid damage and ensure proper functioning.
- 17.** Follow the instructions carefully and wear the goggles correctly to avoid strong light stimulation.
- 18.** Set the irradiation intensity according to your skin condition. Always start with Gear 1 for the first use.
- 19.** It is normal if the skin appears slightly red or warm after use.
- 20.** The device should not be left unattended while turned on. Always turn it off after use.
- 21.** Clean the skin before use to ensure it is free of hair residue, completely dry, and free of oily substances or stains.
- 22.** If there is a significant temperature difference between the operating environment (more than 20 degrees), turn the device off and wait for about 2 hours before using it again.
- 23.** Do not use the device near flammable or explosive objects (such as sprayers, gas stations, etc.) to avoid risks of fire or electric shock injury.
- 24.** The recommended ambient temperature for the device is between 0-45°C. Please ensure the room temperature is within this range before use.
- 25.** Avoid direct ultraviolet radiation on the treated area for 24 hours after treatment.
- 26.** If the built-in fan does not run (i.e., you do not hear the buzzing sound) after turning on the device, do not use it to prevent failures such as smoke or fire.
- 27.** Do not use the device in the bath, shower, or other places close to water. Keep the power adapter dry at all times.
- 28.** It is not recommended to use the product for extended periods. After using it for 1 hour, turn off the device and let it rest for 10 minutes before using it again.
- 29.** If the casing, lens, or power supply is damaged, stop using the product immediately. Do not attempt to repair the product without authorization. It must be repaired by professional personnel from the manufacturer's maintenance department or a similar service.

Restrictions on Use

1. The product should not be used for treating diseases.
2. It is not recommended for use on white or gray hair.
3. Do not use the product on naturally dark skin.
4. Avoid sun exposure before using the device for hair removal.
5. Ensure the lamp window is close to the skin when performing hair removal.
6. Do not use the device in a dark environment.
7. For proper functioning, avoid using the device in direct sunlight.
8. Do not hit or shake the product violently.
9. Always turn off the device before unplugging; failure to do so may damage the product.

Skin Sensitivity Test

Before using the device on a new body area, we recommend conducting a patch test on a small area of skin. Wait for 24 hours after the test to ensure your skin is suitable for IPL treatment.

Light Intensity Settings

Always begin with the lowest light intensity level and choose from the nine available intensity settings based on your comfort level. Reduce the intensity if you experience any pain or discomfort. Individuals with darker skin are at higher risk of skin irritation due to increased light absorption, so lower intensities are recommended. People with dark brown to black skin should not use IPL treatment at all. Please refer to the User Reference Chart below:

Skin color Hair color	white	Like white	creamy white	Light brown	Beige	Dark brown black
Black	✓	✓	✓	✓	✓	✗
Dark brown	✓	✓	✓	✓	✓	✗
Brown	✓	✓	✓	✓	✓	✗
Golden	✓	✓	✓	✓	✓	✗
Gules	✓	✓	✓	✓	✓	✗
white	✓	✓	✓	✓	✓	✗
Grey	✓	✓	✓	✓	✓	✗



Applicable




Not Applicable


Pre-Use Preparation

1. Clean the skin thoroughly to ensure it is hairless, completely dry, and free from oily residues, including perfumes, lotions, deodorants, sunscreen, or other similar substances.
2. Inspect the flash window to ensure it is free from dirt, hair, or debris. Check for any damage or abnormalities. If any dirt is present, clean the window to prevent product malfunction, deterioration, or skin issues.
3. Perform a skin sensitivity test and wait 24 hours to assess the skin condition. If no abnormal reactions, such as redness or blisters, are observed in the treated area, the product can be used safely. Hair removal should not be performed until the skin sensitivity test has been completed.


Instructions for Use

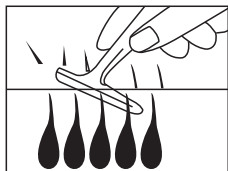


1. Use a razor to shave the hair in the areas that will be treated. Do not touch the skin directly with the razor. It is recommended to use shaving foam or shave while taking a bath. Since the razor is sharp, use it carefully to avoid scratching the skin.

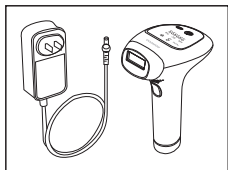


2. Clean the light outlet and connect the power adapter.

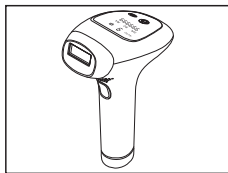


3. Long press or short press the power button to turn on the device. Once powered on, press the ice button to activate the ice cooling mode.

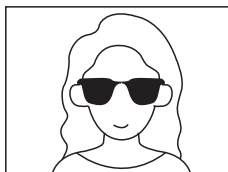

4. Wear goggles to protect your eyes before starting the hair removal operation.

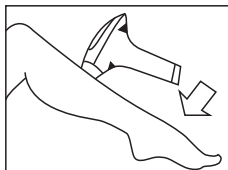

5. Place the device vertically against your skin, press the button once, then move the lamp to the next area without overlapping. Continue operating the device in the treatment area in order.

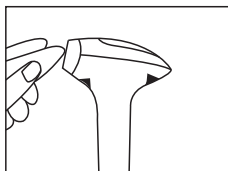

6. It is normal for water droplets to appear on the cooling area of the lamp during the Ice Cooling mode. Please use a soft cloth to wipe off the water droplets after turning off the device.

Hair Removal Tips

- a. The continuous flashes from the hair removal device will generate heat. To enhance the effectiveness of the ice compress, it is recommended to use the ice lamp head to apply compression to the treated area for approximately 30 seconds after hair removal. Then, proceed to the next area for treatment. This will help calm the skin and improve the ice cooling effect.
- b. Press the lamp head firmly against the skin to ensure full contact with the root of the hair follicle. This will allow the hair follicle to absorb more energy, resulting in a more effective hair removal treatment. With the ice cooling care technology, the device effectively cools the skin, reduces heat, and better protects the skin.

Treatment Schedule Guide

1. Use the device 3 times a week during the first 1-3 weeks. During this period, hair growth slows down, and hair decreases by more than 75 % .
2. Use it 2 times a week during the 4-6 weeks. The hair becomes thinner and less noticeable, and the skin becomes smoother.
3. After 6 weeks, use the device 1-2 times per month, selectively, to achieve permanent hair removal.

Note: Please adjust the treatment frequency based on the thickness and hardness of the hair. The gear and frequency can be increased or decreased accordingly. The specific treatment schedule may vary from person to person.

Post-Treatment Care

1. It is recommended to avoid contact with water for at least 6 hours after the hair removal treatment. After this period, apply moisturizing water to hydrate the treated area.
2. Do not apply perfume, masks, or cosmetics containing scrub beads after hair removal.

- 3. Avoid UV exposure for 24 hours after the procedure and apply sunscreen if necessary.
- 4. Artificial tanning is not allowed during hair removal treatments.
- 5. If you require plastic surgery in the treated area after hair removal, please consult your doctor.
- 6. If you experience discomfort, stop using the product immediately and consult a doctor.

Cleaning & Maintenance

- 1. To ensure optimal performance and the longest service life of the product, clean it before and after each hair removal session to prevent long-term accumulation of hair and dust. Failure to clean the product properly may affect its performance.
 - 2. Before cleaning, ensure the device is turned off and unplugged to allow it to cool. Use a soft, clean, and dry cloth to wipe and clean the surface of the light outlet.
 - 3. Store the device in a dry, dust-free environment with a temperature range of 10°C-60°C.
 - 4. Do not wash the device or its accessories in faucets or dishwashers. Avoid using steel balls, abrasive cleaners (toothpaste, etc.), or corrosive liquids (such as gasoline, engine oil, acetone, etc.) to clean the product.
- Note:** If water alone does not clean effectively, you may use a few drops of high-concentration alcohol (no more than 95 %) to clean the parts.

FAQ

Problem	Possible Reason	Solution
Can't turn on	The device is not properly connected to the power source.	Ensure the power is properly plugged in.

Unable to flash	The temperature of the lamp or device may be too high. The lamp head may not be making full contact with the skin.	Allow the device to cool down before use. Place the device on the skin area and wait for the RA icon to display on the screen.
The display does not show	The device may be defective.	Check if the power supply is working correctly. If not, repair or replace the product immediately.
Smelly when flashing	The filter may be dirty or there could be residual hair on the skin surface.	Clean the filter and remove any hair. Reduce the energy level to avoid overheating.
Unbearable pain during hair removal	The area to be treated was not shaved properly.	Before using the device, shave the area to be treated. If shaving causes skin irritation or allergies, trim the hair as short as possible and use the device at a comfortable intensity setting.
	The area to be treated may not be suitable for this product.	Avoid using the device on sensitive areas such as the inner labia, vagina, anus, nipples, areolas, and the face or scrotum.
	The light intensity is set too high.	Lower the light intensity to a more comfortable level.
	The light transmission window is damaged.	If the device is not working correctly, repair or exchange the product immediately.
After hair removal, the area where the hair was removed turns red	Slight redness is normal, harmless, and will subside quickly.	Soothe the skin with aloe vera gel or an ice pack to calm any irritation.

Unsatisfied with the hair removal effect	The light intensity setting may be too low for your needs.	Increase the intensity setting if needed for better results.
	Flashing was not applied to areas close to previously treated areas.	Ensure that the flashes overlap in areas close to each other during treatment.
	This device may not be effective on your hair color.	If your hair is blonde, gray, red, or white, the device may not be effective.
	The device has not been used consistently or sufficiently.	To ensure proper hair removal, follow the recommended treatment schedule.
	You may be reacting slowly to light exposure.	Use the device continuously for up to 5 months as hair regrowth may still occur.
The product becomes hot and humid during use	Prolonged use of the device may cause overheating.	After use, shut down the device and allow it to rest. Ensure proper heat dissipation before the next use.

Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **<https://www.porodo.net/warranty>** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
<https://www.porodo.net/warranty>

Contact Us

If you have any questions about this Privacy Policy, please contact us at:
info@porodo.net

Website: **<https://www.porodo.net/>**

Service Support: **support@porodo.net**

Instagram: **[porodo](#)**