

porodo | CHOOSE THE BEST
Rechargeable Eye Massager
With Heat & Bluetooth Music
SKU:PD-INS5HD



Safety Precautions

Please read and understand all instructions before use and follow them during use. Do not use the product if you have severe eye symptoms, have had eye surgeries, or suffer from serious eye conditions such as glaucoma, cataract, and retinal detachment.

-Clean the product with alcohol-soaked cotton. Do not use water.

-Do not disassemble the product yourself.

-Be mindful of the product's surface, which may become hot during use, especially if you are insensitive to heat.

-This product is not suitable for those who are physically weak, have slow response times, or have mental disorders, unless they are under the guidance or assistance of their guardians. The product is not designed for children and individuals with disabilities to use alone. Keep the product out of children's reach.

-This product is a massager intended for personal health care and cannot replace the effects of medicines and medical devices.

-Do not disassemble, crush, or squeeze the lithium battery in the product, and do not expose it to fire. Stop using the product if the battery swells significantly.

-Do not use the battery if it has been soaked in water.

-Avoid placing the product in environments with high temperatures, as it contains a lithium battery.

-If you have sensitive skin, facial allergic dermatitis, or a reaction to plastic materials, consult a doctor before using the product on your face.

-Remove contact lenses, hairpins, and hair accessories before use.

-Avoid using the product in hot and humid places.

-Avoid exposure to strong sunlight and gases that can cause corrosion.


-Do not use the product while walking or driving, as it obstructs vision.

-Charge the product with a mobile phone charger rated at 5V/1A.


-Keep your eyes closed while using the product.

-Avoid using sharp objects that could scratch or pierce the product.

Product Structure

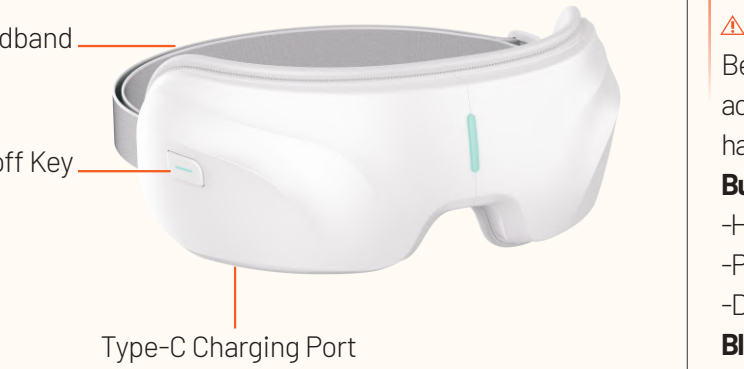

Product Characteristics

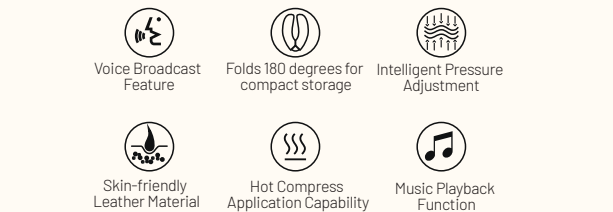

Instructions

Wearing the Product

Position the product over your eyes, secure it, and adjust the strap for optimal comfort. Move the main unit up/down and left/right to find the most comfortable position.

Attention

Before use, remove glasses, contact lenses, hairpins, and any hair accessories to prevent discomfort. Do not use the product with wet hands. Avoid dropping or stepping on the product.

Button Control

-Hold down the power button for 2 seconds to turn the device on or off.

-Press the power button to cycle through modes.

-Double-press the power button to enable mute.

Bluetooth Music

Bluetooth is active by default, playing the device's built-in music upon powering on and connecting to Bluetooth devices. It reverts to built-in music when Bluetooth is disconnected.

(Bluetooth Name: Por-Eyecare)

Note

Extended listening at high volumes may cause hearing loss.

Before Use

1. Fully charge the battery using a 5V mobile phone charger.
2. Use a U8 plug to connect the mobile phone charger with the USB cable.
3. Attach a micro USB plug from the main device to the mobile phone charger to initiate the charge.
4. The battery indicator will flash during charging and remain lit when charging is complete.
5. The device will not operate while charging and will shut down if charging begins during use.

Modes of Use

- *Default: Eye Protection Mode: Air Pressure + Hot Compress + Music(ten minutes).
- *One-touch Loop Mode: Voice broadcasts the current mode after pressing the mode button.
- *Duration: Each mode is active for ten minutes by default, with automatic shutdown afterward.
- *Eye Protection: Air Pressure + Hot Compress + Music
- *Eye Beauty: Air Pressure + Hot Compress + Music
- *Soothing: Air Pressure + Music
- *Nice Sleep: Air Pressure
- *Hot Compress: Hot Compress + Music

Technical Parameters

Product Name	Rechargeable Eye Massager with Heat & Bluetooth Music	Power Supply	Built-in Battery
Control Mode	One-key Operation	Battery Capacity	1200mAh
Message Mode	Five Kinds	Function	Airbag / Hot Compress / Music
Time of Massage	10 Minutes per Session (automatic shutdown)	Charging Mode	Type-C
Power Rating	5W Less than	Charging Current	1A/5V
Product Size	112.5mm*90*221	Net Weight	297g

Hazardous Substances in Products

Part Name	Lead (Pb)	Mercury (Hg)	Cadmium (Cd)	Hexavalent Chromium (Cr(VI))	Polybrominated Biphenyl(PBB)	Polybrominated Diphenyl Ether PBDE
Plastic Parts	○	○	○	○	○	○
Circuit Board	×	○	○	○	○	○
Pressure Control Component	○	○	○	○	○	○
Battery Pack	×	○	○	○	○	○
Straps	○	○	○	○	○	○
Attachments	○	○	○	○	○	○

0: The content of the hazardous substance in all homogeneous materials of the component is below the required limit according to GB/T 26572.

X: The content of the hazardous substance in at least one homogeneous material of the component exceeds the required limit according to GB/T 26572.

Cleaning and Storage

- *Turn off the power before cleaning. Clean the surface of the product with a dry soft cloth.
- *When not in use for extended periods, cover the product with a dust shield(not provided).


Attention

Before cleaning, disconnect the power supply. Avoid using corrosive cleaners like benzene or diluent, which may cause discoloration or damage to the product structure. Do not expose the product to sunlight or high temperatures to prevent damage.

Common Problems and Solutions

	Problems	Possible Causes	Suggested Solutions
1	Noises like "whew" and "da".	Normal working noises from the pump inflating/deflating, motor spinning, the clothing cover rubbing, etc.	Unnecessary.
2	Automatic shutdown after startup.	The lithium battery is low.	Charge the battery.
3	The power button doesn't work.	The battery is not charged.	The charger maybe broken. Charge the battery.
4	Air pressure is too low.	The lithium battery is low.	Charge the battery.

Battery Removal Diagram

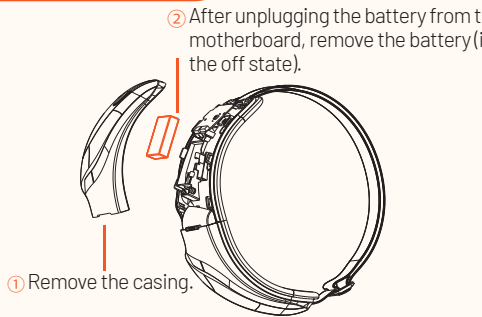

*If issues persist after these solutions, turn off the product and contact the dealer or the company.

*The company is not liable for any loss or damage due to unauthorized repair or dismantling.

Charging Services

- *Damage from improper use or storage by the user is not covered.
- *Damage from service outlets not appointed by the company is not covered.
- *Products outside the warranty period are not covered.

Disposal

1. Remove the battery before discarding the product.
2. Cut off the power supply before removing the battery.
3. Do not replace the battery without authorization.
4. Use the appropriate tools to remove the upper cover of the product, then remove the battery and dispose of it at the appropriate recycling facilities.

Transportation and Storage

- *Ambient temperature: -40°C to 50°C
- *Relative humidity: less than %80 with good ventilation. Avoid corrosive gases.
- The product is suitable for transport by air, rail, road, and sea. It should not be exposed to rain or snow, nor should it be stored upside down or crushed.
- If the product has been stored for over three months, inspect the product and power it on for half an hour before use.

Warranty

Exclusion

- Porodo doesn't cover any damage caused by the normal use of the product over time during the warranty period.
- The Porodo warranty doesn't include coverage for improper use, deliberate damage, accidents, theft, carelessness, or any repairs done by unauthorized service centers.
- The power needs for each product are clearly stated on their packages. Please check these requirements against those of your country before buying. Porodo won't be responsible for any problems related to this.

For more info, please check:
porodo.net/warranty

Contact Us

If you have any questions about this Privacy Policy, please contact us at info@porodo.net

Website: porodo.net

Service Support: support@porodo.net

Instagram: porodo