

PORODO
LIFESTYLE



Porodo Lifestyle
Handheld Hot-Stone Massager

SKU: PD-LSAMRC

Table of Contents

Features	2
Specifications	2
Usage Prohibitions	2
Schematic View	3
Instructions for Use	3
Charging	3
Warnings	4
Battery Disposal	4
Warranty	5
Contact Us	5

Features

1. The 10 intensity levels allows users to customize the intensity of the massage to suit personal comfort and therapeutic needs.
2. Utilizes an acupuncture stone to enhance massage strength, deepen muscle relaxation, and improve relief from discomfort.
3. Integrates heat for a warming compress, enhancing the relaxation and therapeutic benefits of the massage
4. The Low Frequency Vibration provides gentle vibrations to soothe and relax muscles.
5. The Medium Frequency Vibration delivers moderate vibrations for targeted muscle relief.
6. The High Frequency Vibration offers intense vibrations to effectively reduce muscle fatigue and penetrate deep into muscle tissues for quick relaxation.

Specifications


Product Name	Porodo Lifestyle Handheld Hot-Stone Massager
SKU	PD-LSAMRC
Charging Voltage	5V
Rated Voltage	7.4V
Charging Current	2A
Rated Power	24W
Battery Capacity	1800mAh
Size	8.5×8×19cm
Net Weight	0.54kg

Usage Prohibitions

1. Do not use the device if you have any of the following conditions:
 - 1.1 Aneurysms, active bleeding disorders, or if you are taking blood thinners.
 - 1.2 Heart disease, or if you have a pacemaker or defibrillator implanted.
 - 1.3 Current pregnancy or diagnosed cancer.


2. Avoid using the device within a 3-inch radius of any surgical site for 90 days post-orthopedic surgery.
3. Do not use the device on sensitive areas including the head, face, cervical spine, or spinal vertebrae.
4. Avoid using the device on areas with minimal muscle coverage near bones, such as the tibia, dorsum of feet, and back of hands.

Schematic View

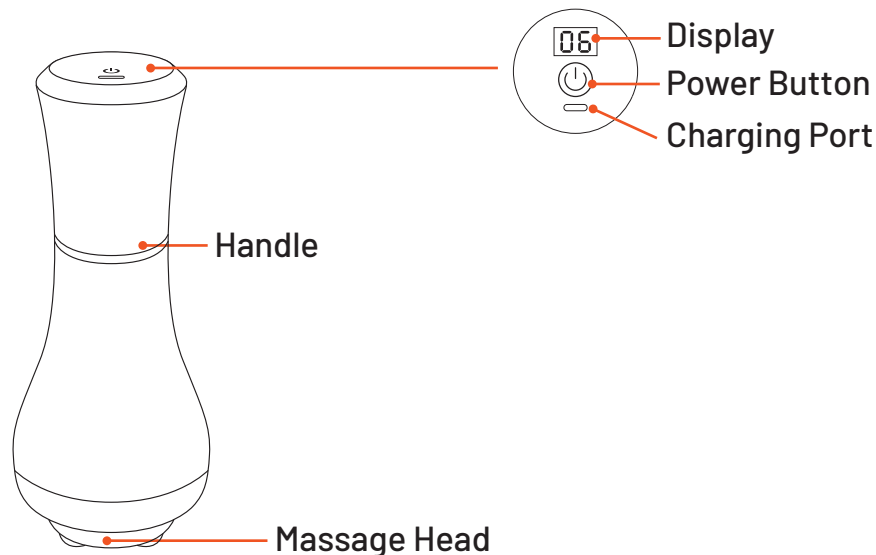


Instructions for Use

1. The device starts in the first gear.
2. Press the power button once to increase the setting by one level.
3. The device has ten intensity levels. Upon reaching the tenth gear, pressing the power button again will reset it to the first gear. The current gear will be displayed digitally.

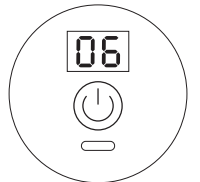

Charging

1. First Charging: Prior to the first use, ensure the battery is fully charged.
2. Charging Indicator: The indicator light will flash during charging and will remain lit once the battery is fully charged.

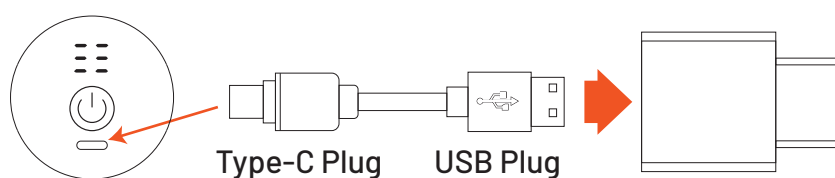

Warnings

- 1.** Do not drop or collide with the charger. If the charger (adapter) is damaged, immediately contact an authorized service center for a replacement.
- 2.** Avoid touching the power cord with wet hands, and do not disconnect the charger (adapter) by pulling on the power cable to prevent electric shock.
- 3.** If the charger is exposed to rain, liquids, or excessive moisture, cease using it and consult an authorized service center for replacement.
- 4.** Keep the battery away from high temperatures and nearby heating devices like sunlight exposure, heaters, microwaves, ovens, or water heaters to prevent the risk of overheating or explosion.
- 5.** Do not disassemble, modify, or immerse the battery in water or other liquids.
- 6.** Avoid inserting foreign objects into the battery compartment to prevent battery leakage, overheating, or explosion.
- 7.** Remove the battery from the device before disposal.
- 8.** Ensure the appliance is turned off before battery removal.
- 9.** Dispose of batteries safely in accordance with local regulations.
- 10.** Charge the appliance using a 5V=2A (10W) adapter that complies with the GB2005-4706.1 standard.
- 11.** This appliance is not suitable for use by unsupervised children.
- 12.** This product should not be used by individuals with physical disabilities, mental or intellectual disabilities, or those lacking the experience and knowledge required to operate it safely.

Battery Disposal

- 1.** Ensure that the product is powered off before proceeding with battery removal.
- 2.** Authorized personnel must use a screwdriver to open the casing of the product. Once the casing is removed, the battery can then be safely extracted.
- 3.** After the battery has been removed, it should be turned over to a designated battery recycling department or deposited in a battery recycling box for environmentally responsible disposal.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **porodo.net/warranty** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
porodo.net/warranty

Contact Us

If you have any questions about this Privacy Policy, please contact us at:
info@porodo.net

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **porodo**