

PORODO
LIFESTYLE



Porodo Lifestyle
W-SONIC Electric
Toothbrush

SKU: PDLFST200BK

Table of Contents

Important safety information	2
Important Usage Prohibitions	2
Schematic View	3
Specifications	3
Instructions for Use	4
Function and Mode	4
Charging the Product	5
Danger Warning	6
Maintenance	6
Environmental Protection Statement	7
Warranty	7
Contact Us	7

Important safety information

1. Please read this important information carefully before using the product and keep it for future reference.
2. If you have had oral or gum surgery within the past two months, consult your dentist before use.
3. If you have a heart rate modulator or any other medical device, consult your doctor or the device manufacturer before use.
4. Only use toothbrush heads recommended by our company.
5. This product is designed solely for cleaning teeth and gums and should not be used for any other purpose.
6. This product is not suitable for children, infants, or anyone unable to operate an electric toothbrush.
7. Keep the product out of reach of children and babies, and do not use it as a toy.
8. To prevent bacterial growth on the bristles and maintain oral health, replace the brush head every three months.
9. When the bristles appear worn or frayed, it is recommended to replace the brush head promptly.


Important Usage Prohibitions

1. Do not use the product if it is obviously damaged before use.
2. This product contains a built-in rechargeable battery. Do not put it into fire, heat it, charge it improperly, or use it at high temperatures.
3. It is forbidden to break, modify, disassemble, or repair this product by yourself.
4. Charge this product only within the temperature range of 0°C to 40°C.
5. When removing the battery, make sure the product is powered off.
6. Batteries should be handled with care and disposed of properly.

Note:

It is normal to experience slight tingling or minor bleeding in healthy teeth during the first use of the electric toothbrush because your teeth are not yet accustomed to this stimulation. This condition typically stops within 2–1 weeks. If it persists beyond 2 weeks, please discontinue use and consult your dentist.

Schematic View

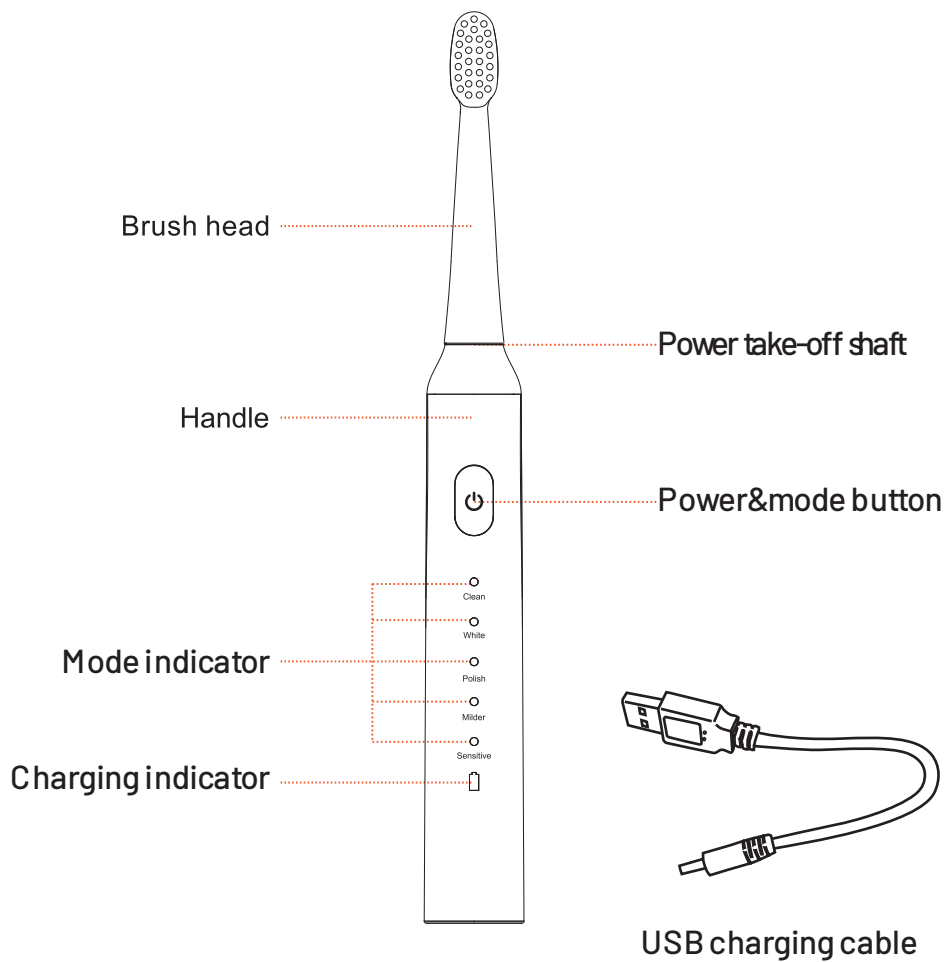

Specifications

Input	5V/1A
Battery Capacity	500mAh
Battery Life	Up to 50 Days
Charging Time	3H
Noise Level	<55dB
Water-Resistance Level	IPX7
Charging Type	DC Port
Vibration Frequency	18000-20000 Stroke/Minute
Product Size	24 x 25x 235mm
Age Group	Adult
Model Number	LFS200

Instructions for Use

1. Insert the brush head vertically into the gap between the brush head and the handle.

Note: A small gap will remain between the brush head and the handle.





2. Wet the brush head and apply an appropriate amount of toothpaste.

3. Gently press the brush head against the junction of your teeth and gums at a 45-degree angle.

4. Press the mode button to select the desired brushing mode. Gently move the toothbrush head back and forth between your teeth.

5. The toothbrush operates intelligently for 2 minutes, pausing every 30 seconds to remind you to move to the next area.

6. After brushing, rinse the toothbrush head with clear running water and shake off any excess water.

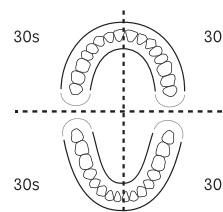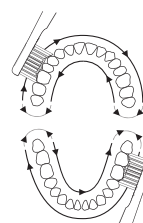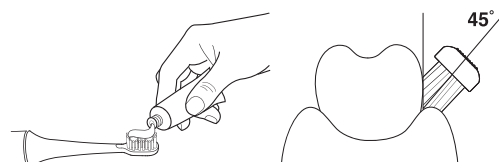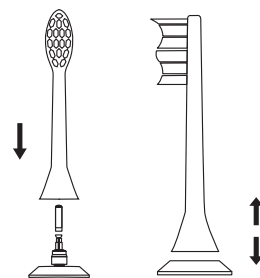

Function and Mode

Short press the power button to turn on the device, long press to turn it off. After turning on, short press to switch between cleaning modes. When the last mode is reached, short press to turn off.

Clean Mode

Effectively removes dental plaque (default mode). The clean mode lasts 2 minutes and includes four 30-second intervals.

White Mode

Whitens teeth in 2 minutes by effectively improving the color of the tooth surface.

Polish Mode


Brightens and polishes the front teeth. Users can press and hold the power button to stop polishing early according to their needs.

Milder Mode

Recommended for beginners. This mode effectively cleans teeth and gently massages gums to improve gum health.

Charging the Product

1. Direct charging: Insert the DC connector of the charging cable into the device's charging port, and connect the other end to a standard 5V $\overline{\sim}$ 1A adapter. (Refer to Figure 1)
2. When the battery power is low, the charging indicator light will flash, indicating that the device needs to be charged.
3. During charging, the indicator light will flash slowly and continuously.
4. Once fully charged, the indicator light will stay on steadily and stop flashing.

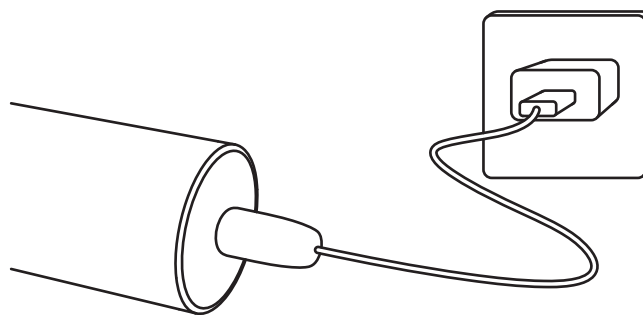

(Figure 1)

When charging, it is normal for the handle and base to feel slightly warm. To ensure safety, always use the original standard charger.

After charging, please unplug the charger to prevent overheating, which can cause component aging, damage the charger, or pose an electric shock risk.

Before connecting the charger to a power source, verify that the voltage on the product matches your local power supply voltage.

Danger


Keep the charger away from any water sources. Do not place or store it near containers holding water, such as bathtubs, washbasins, pools, etc. Never immerse the charger in water or any other liquids. After cleaning, ensure the charger is completely dry before plugging it into the power supply.

Danger Warning


1. Do not replace the power cord. If the power cord is damaged, discard the charger immediately. Always use the original charger model to avoid safety hazards. If any part of the product (brush head, toothbrush handle, or charger) is damaged, stop using it immediately.
2. Use an adapter with an output voltage of 5V and 1A to charge the product. Charging with an adapter that has excessive output voltage may cause damage or pose safety risks.

Maintenance

1. When cleaning the brush head, please turn it downward.
2. After each use, clean and air-dry the bristles and brush heads naturally.
3. Avoid pulling the bristles to prevent deformation.
4. When cleaning the toothbrush body, wipe it with a dry cloth.
5. Do not soak the toothbrush in water for extended periods.
6. Avoid using water or acidic cleaning solutions above 40°C for cleaning.

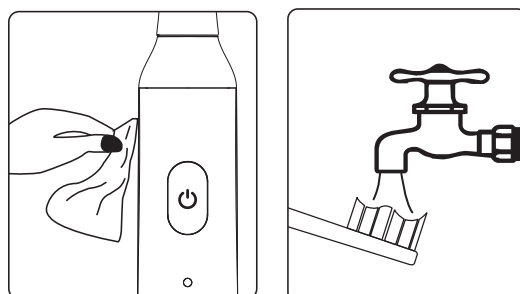

Environmental Protection Statement

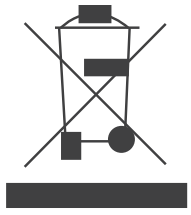

This sonic electric toothbrush contains a built-in rechargeable lithium-ion battery. When disposing of this product, please do not discard it with regular household waste. Remove the battery and take it to a government-designated battery recycling center to help protect the environment.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **porodo.net/warranty** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
porodo.net/warranty

Contact Us

If you have any questions about this Privacy Policy, please contact us at:

info@porodo.net

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **[porodo](https://www.instagram.com/porodo)**