

PORODO
LIFESTYLE



Porodo Lifestyle

**ULTIMATE ORAL CARE DUO Electric
Toothbrush and Water Flosser Combo**

SKU: PDLSTC2254

Table of Contents

Preface	2
Warning	2
Water Flosser Components	3
Water Flosser Specifications	4
Instructions for Using the Water Flosser	4
Product Charging Guide	6
Care and Maintenance	8
Toothbrush Components	9
Toothbrush Modes	9
Instructions for Using the Toothbrush	10
Toothbrush Charge Indicator	11
Toothbrush Specifications	12
Disposal	12
Warranty	13
Contact Us	13

Preface


Thank you for choosing The Ultimate Oral Care Duo Electric Brush and Water Flosser for the perfect smile. We trust this combo will meet all your dental needs and deliver the desired results. Please review the manual before using the product to ensure safety and optimal performance.

Warning

- 1.** Do not dismantle the product or alter any components without proper authorization, as it may cause damage or electrical shock, invalidating the warranty.
- 2.** Avoid cleaning the product with water above 40°C. If the product is submerged during charging, unplug the power adapter and remove the product immediately to prevent electric shock.
- 3.** Do not expose the product to fire or excessive heat. Avoid charging the product in hot or humid environments and ensure it is not placed in direct sunlight for prolonged periods.
- 4.** Do not unplug the power plug or charging cable with wet hands while charging.
- 5.** The product does not include an adapter. If using a third-party adapter, ensure it outputs DC5V, similar to a mobile phone charger, and verify it meets the required voltage specifications.
- 6.** If the product will not be in use for an extended period, unplug both the power plug and charging cable to avoid potential hazards.
- 7.** This product is not suitable for use by children or individuals with limited physical, sensory, or mental capabilities, unless supervised by someone with the relevant knowledge to ensure safe operation.
- 8.** Avoid using the product if you have oral ulcers, recent oral surgery, sensitive gums, or teeth, unless approved by a dentist. Use may cause bleeding or additional damage.
- 9.** If you experience any bleeding while using the product, consult your dentist.
- 10.** Always ensure the nozzle is properly inserted before activating the water flosser to prevent unnecessary splashing.

- 11.** Dispose of Li-ion batteries responsibly and in accordance with local electronic waste disposal regulations.
- 12.** Do not charge the product directly or with the charging base when wet or in the shower.
- 13.** Avoid exposing the product to direct sunlight.
- 14.** The product is waterproof and can be used while bathing or showering.
- 15.** If you have a pacemaker, consult your doctor before use to avoid interference.
- 16.** Ensure the charging port is dry and free from any water before charging the product.

Water Flosser Components

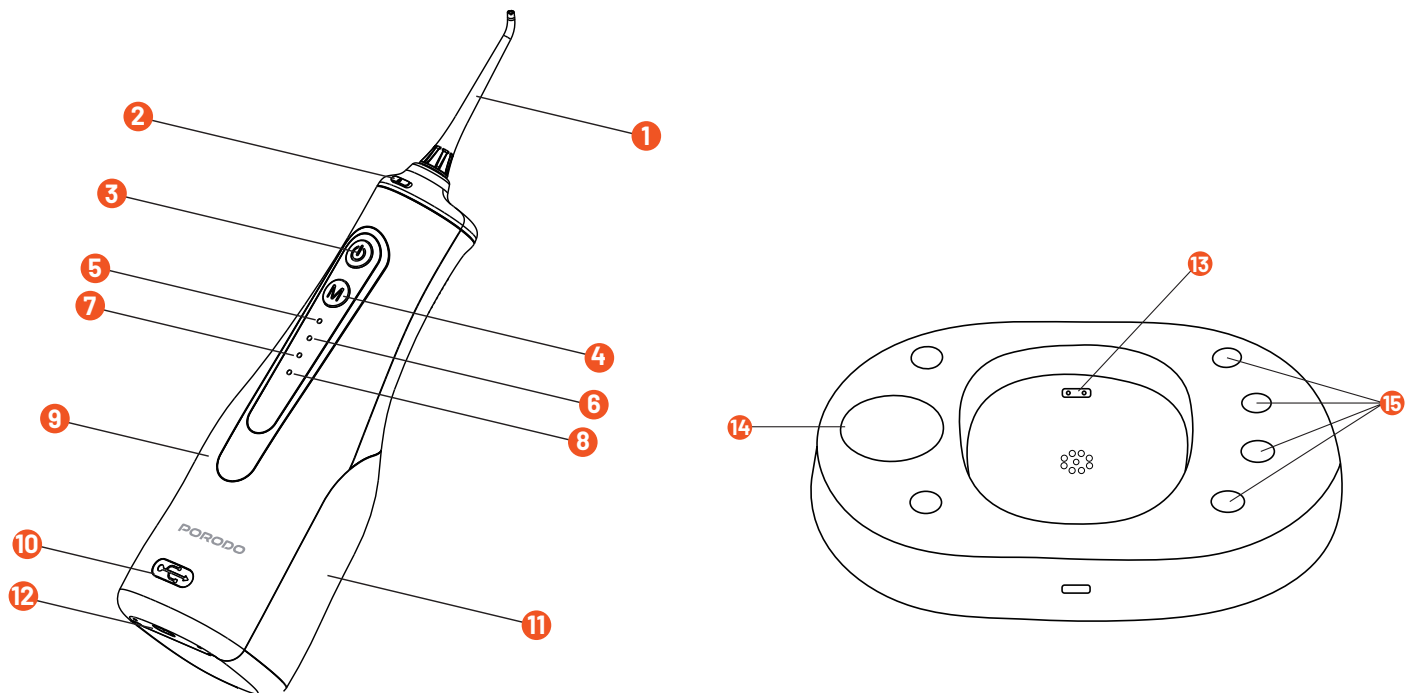

- | | |
|---------------------------------|---|
| 1. Nozzle | 9. Main body |
| 2. Nozzle release key | 10. Type-C charging port |
| 3. Power button | 11. Water tank |
| 4. Mode button | 12. Charging port |
| 5. Pulse mode indicator | 13. Water flosser charging base and holder |
| 6. Soft mode indicator | 14. Sonic electric toothbrush holder |
| 7. Normal mode indicator | 15. Flosser nozzles holder |
| 8. Power mode indicator | |


Water Flosser Specifications

Rated Power	6W
Battery Capacity	1500mAh
Charging Time	4H
Working Time	2 Weeks (2x/Day)
Charging Interface	Type-C & Charging Base
Water Tank Capacity	300mL
Product Size	69 x 85×224mm
Flosser Modes	4 Modes (Pulse/Soft/Normal/Power)

Instructions for Using the Water Flosser

1. Nozzle Mounting


Insert the nozzle into the nozzle holder located at the top of the main unit until it clicks into place.

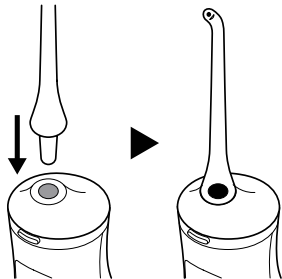

2. Injection of Cleaning Solution

2.1 Open the lid of the water tank filling port, pour clean water into the tank, and then close the water tank lid.

2.2 If using a concentrated mouthwash, dilute it with water before filling the tank.

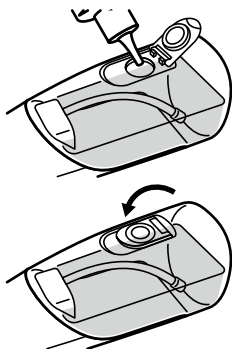

Usable Cleaning Solutions


The following liquids are recommended for use in the water tank

Important: Do not use any other liquids, as they may cause malfunction.

1. Tap water
2. Warm water (40°C or below)
3. Mouthwash


3. Mode Settings

Press the mode switch button on the front of the main unit to select the desired water pressure. Each press cycles through the modes in the following order: Pulse, Soft, Normal, and Power. Check the corresponding mode indicator lamp to confirm the selected mode.

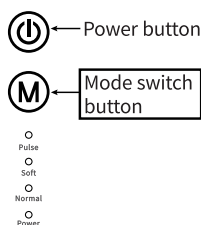

4. Power On

Insert the tip of the nozzle into your mouth and press the power button. The indicator light will illuminate to signal that cleaning has begun.

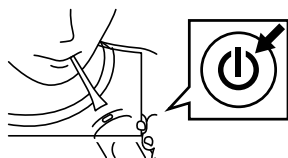


5. Cleaning

5.1 While cleaning, lean forward and slightly open your mouth to allow the cleaning water to be flushed out.

5.2 Please note that if you open your mouth too wide, the cleaning water may splatter.


5.3 The main unit is equipped with a waterproof mechanism, ensuring safe use even if it gets wet during cleaning.

5.4 The nozzle can be rotated 360°, allowing you to adjust it to a comfortable angle for use.

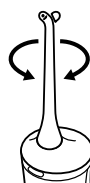


Note: If the water supply hose inside the tank rises above the water surface, the cleaning water will not flow out. Adjust the main body and ensure the hose is submerged to allow proper water flow.

Caution: Aim the water stream between the teeth and gums, then move slowly along the gums from the molars (back teeth) to the front teeth. Clean the surfaces and the back of the teeth.

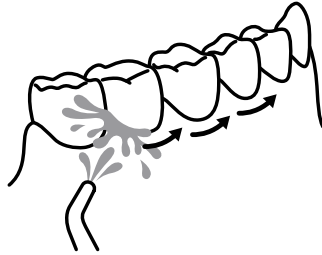


6. Stop

To stop the operation of the main unit, such as when the cleaning water is depleted, press the power button again.

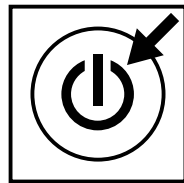

7. Remove the Nozzle

Press the nozzle release button and pull the nozzle upward to remove it.

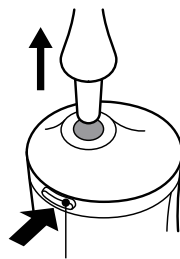

Product Charging Guide

If you stop or press the power button during operation, and all four mode indicators flash simultaneously, the battery power is low and needs to be charged.


1. Remove the cap from the charging port of the main unit and connect the supplied USB-C cable to the port. Then, connect the other end (USB-A) of the charging cable to your USB output charging adapter.

1.1 The water flosser can also be charged when placed on the charging base.

1.2 If the charging port of the main unit is wet with water, wipe it with a dry cloth before connecting the charging cable. Only start charging after the port is completely dry.


2. Charging

When charging begins, the four mode indicator lamps light up sequentially to indicate that charging is in progress.

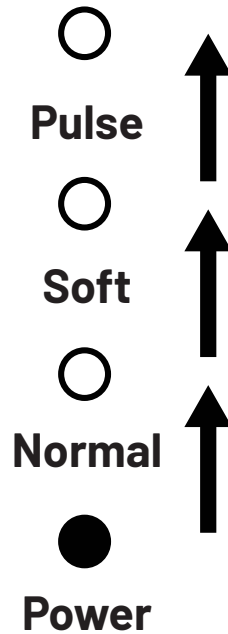


3. Fully Charged

Once charging is complete, all four mode indicator lamps will light up simultaneously, signaling that the device is fully charged.

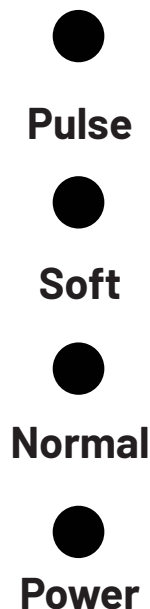


Care and Maintenance

1. After each use, open the tank and discard any remaining water. Then, press the power button to activate the unit and drain any water remaining in the nozzle. When the water stops flowing, press the power button to turn off the device.
2. If dirt accumulates inside the tank, add a small amount of lukewarm water (40°C or lower), swish it around, and then perform the above steps.
3. To maintain cleanliness, rinse the nozzle with water after each use and allow it to dry.

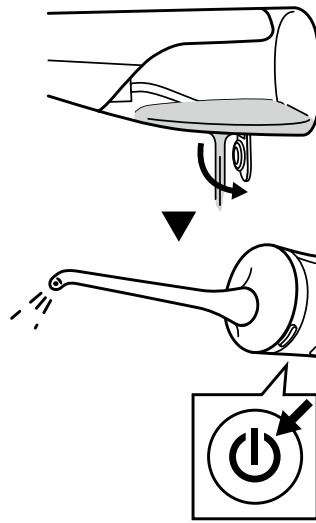


If no water sprays for an extended period, unplug the nozzle and press the water valve in the water outlet intermittently with a tool that has a diameter of less than 6mm until water begins to spray.

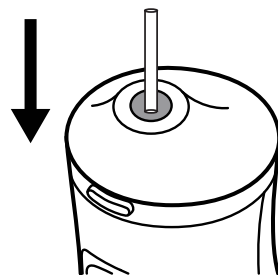


Slide the water tank to detach it from the base.

There is an air inlet in this section of the water tank. When the tank is filled with water and turned upside down, a small amount of water may flow out, which is normal. During startup, if air enters the tank, ensure the air pressure inside is balanced.

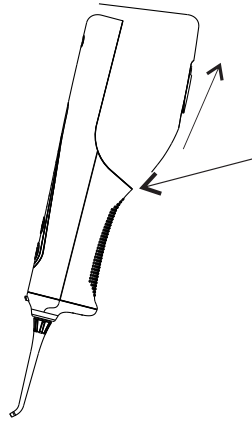

Toothbrush Components

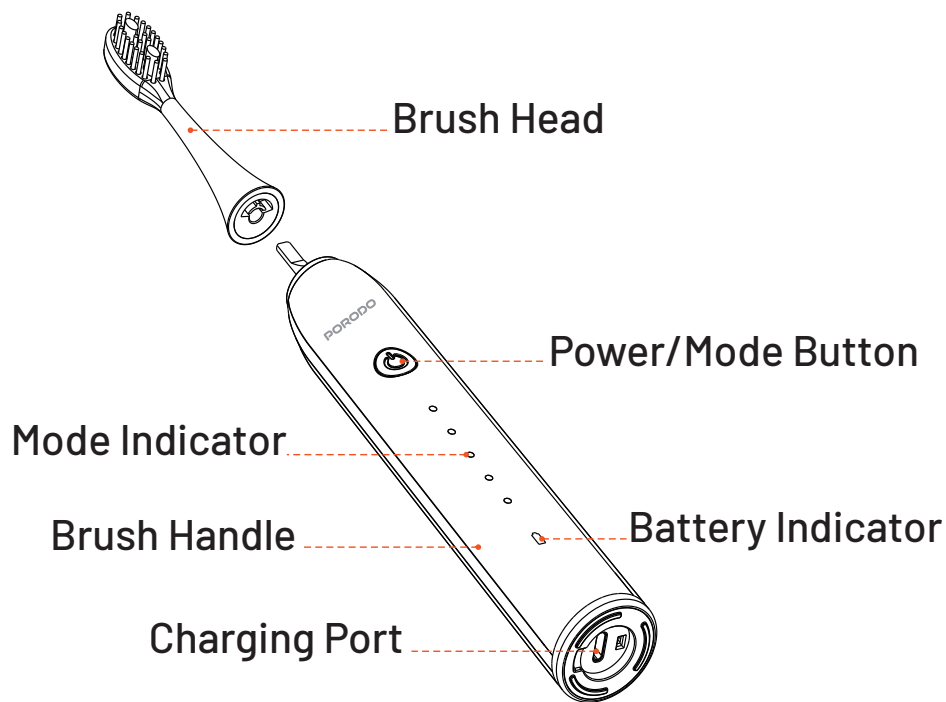


Toothbrush Modes

- Clean:** Removes food particles and dental plaque.
- Refresh:** Rinses the mouth and refreshes the teeth.
- White:** Removes surface stains and whitens the teeth.
- Polish:** Removes calculus and tartar from the teeth.
- Sensitive:** Massages sensitive gums and reduces vibrations for a gentler experience.




Instructions for Using the Toothbrush


1. Attach the brush head firmly to the handle shaft, leaving a gap of approximately 1mm between the brush head and the handle. This gap allows effective vibrations for optimal performance of the toothbrush.

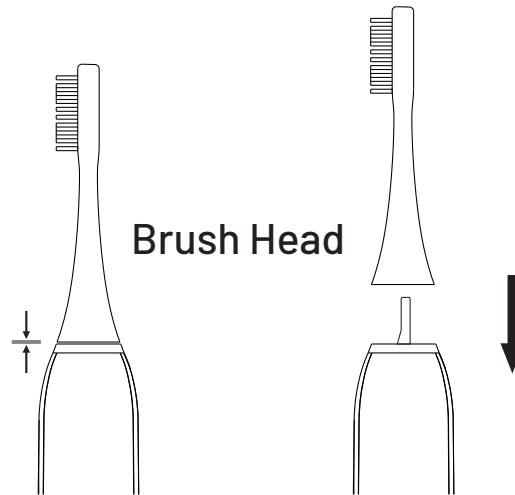


2. Position the brush head at the border of your teeth and gums.

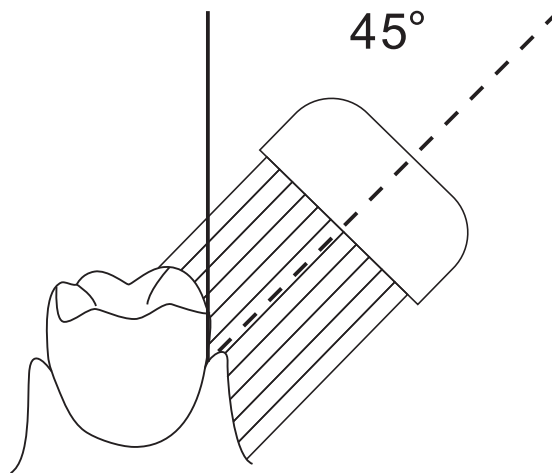


3. Place the brush head gently at the junction of the teeth and gums. Press the mode button several times within three seconds to select the desired operation mode. Move the brush head gently back and forth to clean all of your teeth.

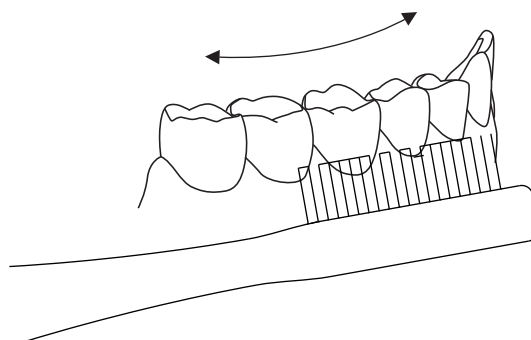

4. Press the power/mode button to select your desired mode. Gently move the brush head back and forth to clean all of your teeth. The toothbrush has a 30-second interval pause to signal when to move to the next section of your mouth. It also features a dentist-recommended 2-minute automatic power-down function.

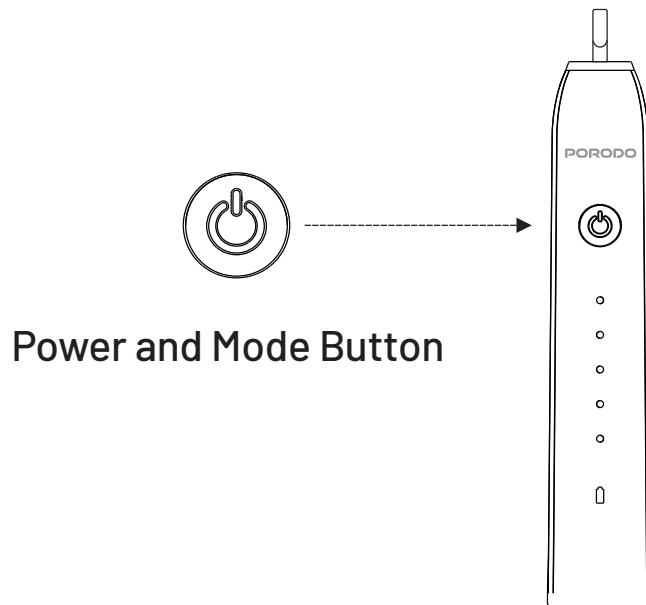


5. Rinse your mouth and brush head with water. Keep the brush head dry and clean the handle with a soft cloth.

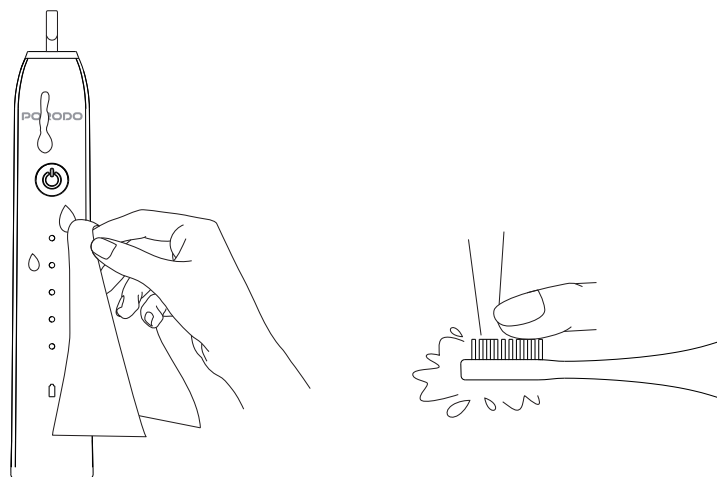

Toothbrush Charge Indicator

Plug the USB charging cable into a USB charger or computer USB port. The toothbrush indicator light will change from red to green after 6 hours, indicating that charging is complete.

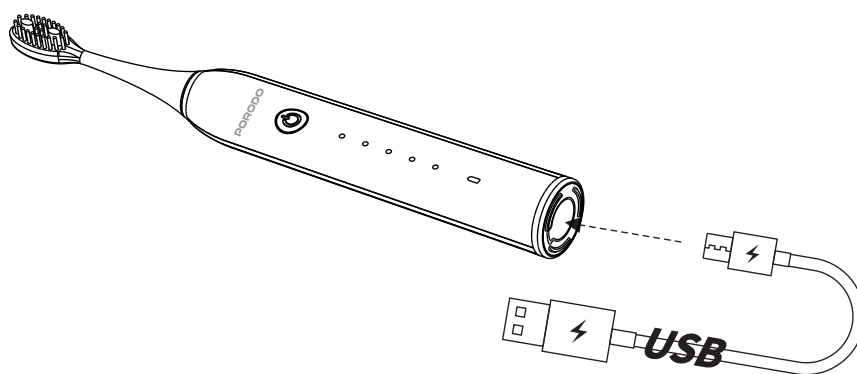

Toothbrush Specifications

Input	5V/1A
Rated Power	1.5W
Battery Capacity	1200mAh
Charging Time	2H
Working Time	2 Months
Charging Interface	Type-C
Noise Level	<60dB
Water-Resistance Level	IPX7
Product Size	32 x 32× 235mm
Toothbrush Modes	5 Modes (Clean/Refresh/White/Polish/Sensitive)

Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **porodo.net/warranty** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
porodo.net/warranty

Contact Us

If you have any questions about this Privacy Policy, please contact us at:
info@porodo.net

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **porodo**