

PORODO
LIFESTYLE



Porodo Lifestyle
W-SONIC Electric Toothbrush

SKU: PDLFST905


Table of Contents

Product Overview	2
Schematic View	2
Specifications	3
Package Contents	3
Charging guideline	3
Battery Indicator Guide	4
Electric Toothbrush Guide	4
Safety Precautions	6
Hazardous Substances in the Product	7
Disposal	7
Warranty	8
Contact Us	8

Product Overview

Our electric toothbrush is engineered to provide a scientifically-driven approach to oral health, combining cutting-edge technology with effective care solutions. Recognized and trusted by countless consumers for its superior performance, this toothbrush delivers a deeper, more thorough clean compared to traditional brushing. It is designed to promote healthier gums, reduce plaque buildup, and enhance overall oral hygiene. With multiple brushing modes and advanced features, it is tailored to meet various dental needs, ensuring a customized experience. To ensure optimal performance and longevity of the product, we highly recommend reading the detailed instructions for proper usage and maintenance. This will help you get the most out of your electric toothbrush, ensuring your oral care routine is as efficient and effective as possible.

Setup Guide

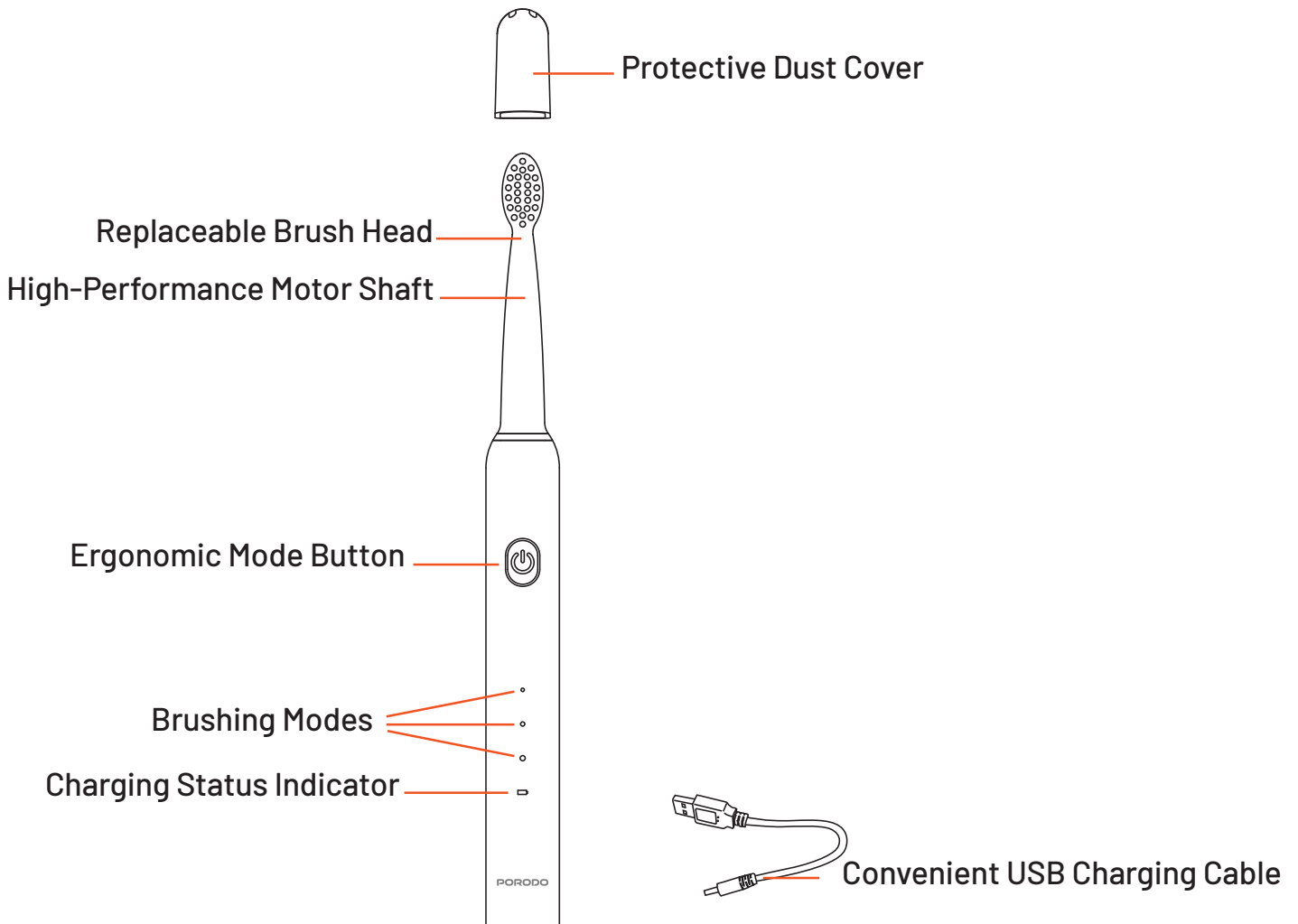

Specifications

Input Power	5V/1A
Battery Capacity	500mAh
Battery Life	Up to 100 Days
Charging Time	3H
Noise Level	<55dB
Water Resistance Rating	IPX7
Charging Port Type	Type-C
Vibration Frequency	19000 Stroke/Minute
Product Dimensions	2.1×22.7cm
Net Weight	51g
Model ID	LFS905

Package Contents

Rechargeable Electric Toothbrush	1
Replacement Brush Heads	2
User Manual	1
USB Charging Cable (Type-C)	1

Charging Guideline

1. Before using the electric toothbrush, verify that the power supply voltage and the power socket are both functioning properly and within the normal operating range. Confirm that the voltage is compatible with the device to prevent any electrical hazards.
2. Before plugging the electric toothbrush into the power supply, ensure that the USB cable is completely dry. Check that the charging plug is free from any water stains, moisture, or physical damage. Only connect the device to the power source if all components are in a safe and proper condition.
3. When using the electric toothbrush for the first time, power it on to check for initial functionality. If the battery level is low or if there is no power, fully charge the

toothbrush before its first use. This will ensure optimal performance and longevity of the device.

Battery Indicator Guide

1. Indicator status during use	Description
1.1 Indicator on	1. Battery Status
1.2 Red indicator on	1.1 Low Battery, Please Charge
1.3 Indicator not on	1.2 Ultra Low Battery (Product may not function properly)
2. Indicator status during charging	2. Charging Status
2.1 Red indicator on	2.1 Charging
2.2 Indicator turns from red to green	2.2 Charging Complete


Notes:

- 1.** Always monitor the battery indicator during use to ensure that the toothbrush has sufficient power for effective operation.
- 2.** Charge the toothbrush when the battery level is low to maintain optimal performance.
- 3.** If the toothbrush has not been used for a prolonged period, it is recommended to fully charge the device before the next use. This ensures the battery is properly calibrated and helps maintain the health of the battery.
- 4.** To ensure the safety of the charging process and enhance the product's lifespan, please follow the charging instructions provided. Proper charging habits contribute to extending the overall service life of the toothbrush.

Electric Toothbrush Guide

- 1.** When the power is low and the LED indicator flashes red, please charge your toothbrush. Refer to the charging instructions for details on how to charge it properly.
- 2.** Attach the brush head to the handle as shown in the diagram. There should be a slight gap between the brush head and the handle, which is completely normal.

3. Apply a small amount of toothpaste to the bristles. Make sure to use a proper amount that will not cause excessive foam.
4. Place the brush head in your mouth, ensuring that the bristles are at a 45° angle to your gums. Turn the toothbrush on by selecting your preferred brushing mode.
5. Move the toothbrush gently up and down, without pressing too hard or scrubbing. Let the toothbrush do the work and avoid vigorous movements for the best results.
6. After turning off the toothbrush, remove it from your mouth to prevent toothpaste foam from splashing. Ensure that you switch the toothbrush off before taking it out.
7. After brushing your teeth, rinse the brush head thoroughly with water to remove any toothpaste residue.
8. Place your toothbrush in a well-ventilated and dry area to ensure proper hygiene and prevent any moisture buildup.

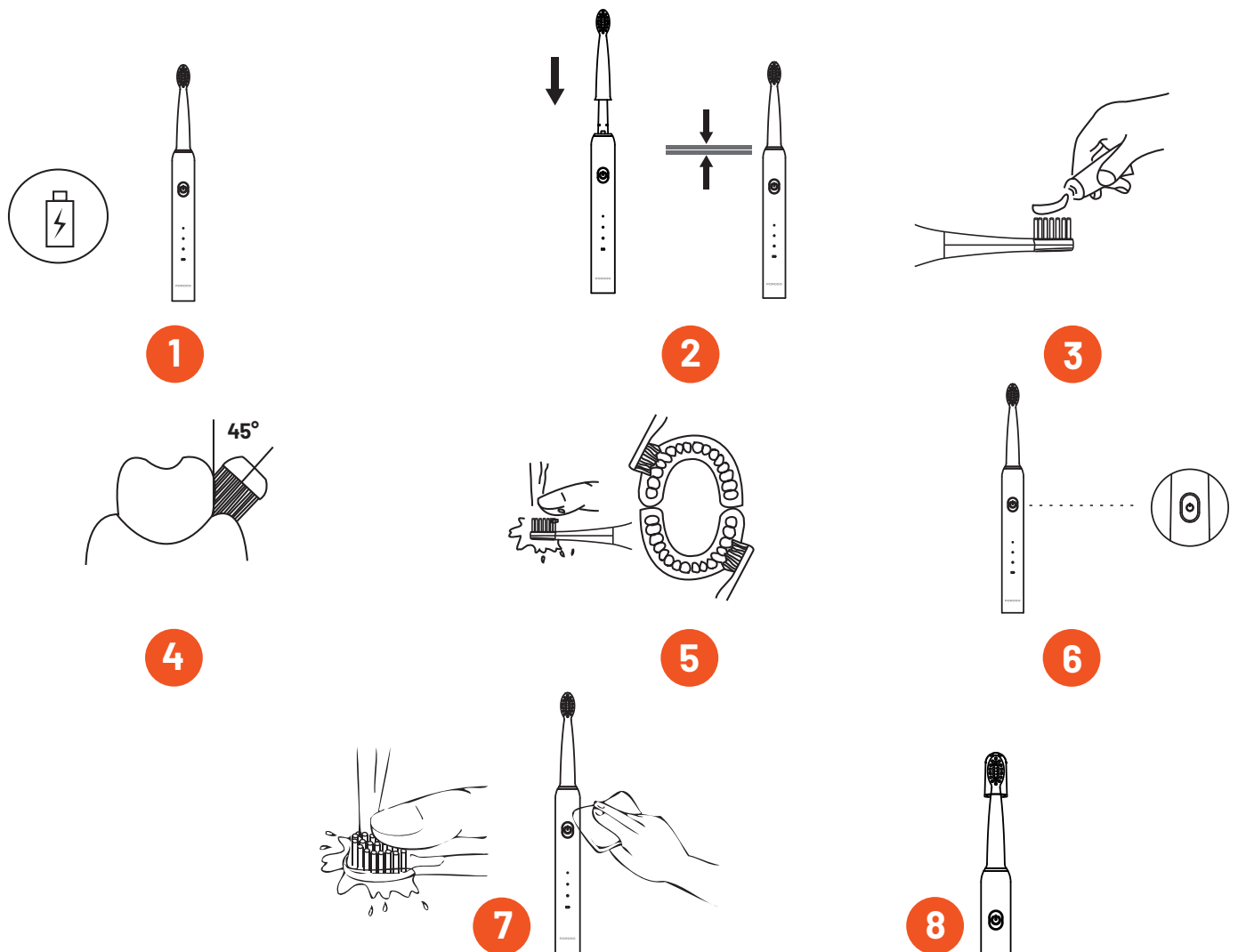

Safety Precautions

- 1.** Only use toothbrush heads that are recommended by the manufacturer.
- 2.** Never expose the product to fire, excessive heat, or high temperatures. Do not charge or use the product in such conditions.
- 3.** For optimal oral health, brush your teeth with your electric toothbrush at least twice a day. Be sure to replace the brush head every 2-3 months to maintain the effectiveness of the brushing action and ensure a thorough clean. Regular replacement of the brush head also helps to prevent the buildup of bacteria, ensuring your brushing experience remains hygienic and efficient.
- 4.** The toothbrush contains a built-in rechargeable battery. Handle with care and follow safety guidelines.
- 5.** Keep the charger away from any water sources. Do not place or store it above or near any container with water, such as a bathtub, washbasin, or pool. Do not immerse the charger in water or any other liquid.
- 6.** Always use the original charger model for charging to avoid potential hazards.
- 7.** Always turn the brush head downward when cleaning it to prevent water from entering the handle.
- 8.** After each use, clean the bristles and brush head, and allow them to dry naturally.
- 9.** Avoid pulling on the bristles, as this may cause them to deform.
- 10.** Dry the toothbrush handle with a dry cloth. Do not soak the toothbrush in water for extended periods. Avoid using water or acidic cleaning liquids above 40°C to clean the toothbrush.

Hazardous Substances in the Product

Components and parts	Hazardous Materials					
	Pb	Hg	Cd	Cr(VI)	PBB	PBDE
Brush Head & Enclosure	○	○	○	○	○	○
PCBA	X	○	○	○	○	○
Battery	X	○	○	○	○	○
Other parts	○	○	○	○	○	○

This form complies with the SJ/T11364-2014 regulations.

(o) Indicates that the content of hazardous substances in all homogeneous materials of the parts is below the GB/T 26572 prescribed limit.

(x) Indicates that the content of hazardous substances in at least one of the homogeneous materials of the parts exceeds the GB/T 26572 limit requirements.

The hazardous substances listed in the table are present in the electric toothbrush but will not harm your health or the environment when used as intended. However, the information about these substances has been left out to help consumers and recycling professionals properly dispose of the product when it is no longer needed. This ensures the product is safely discarded in an environmentally friendly way

Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations. Follow these steps for proper disposal:

- 1.** Ensure the battery is completely drained, so the toothbrush can no longer be turned on. This helps prevent potential electrical hazards during disposal.
- 2.** Use a screwdriver to pry off the bottom cover of the toothbrush, exposing the internal components.
- 3.** Gently press down on the motor shaft and carefully remove the internal parts of the device.

4. Use appropriate tools to sever all connecting wires, ensuring the battery is fully detached. Properly dispose of the battery and wires separately as hazardous materials.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **porodo.net/warranty** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
porodo.net/warranty

Contact Us

If you have any questions about this Privacy Policy, please contact us at:

info@porodo.net

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **porodo**